




neodisher® MultiZym



Detergent for the manual reprocessing of thermostable and thermolabile instruments

Liquid concentrate

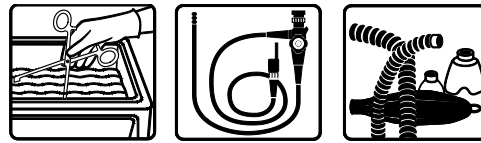

Fields of application:

- Intended for use for the manual cleaning of thermostable and thermolabile instruments, incl. MIS and micro instruments, dental instruments, flexible endoscopes, anesthesia utensils, containers and other medical utensils in immersion and ultrasonic baths

Performance spectrum:

- Reliably removes residues such as dried and denatured blood, fats and secretions
- Causes depletion of organic materials and prevents the redeposition of protein residues
- Excellent material compatibility: suitable for stainless steel, instrument steel, optics, commonly used synthetic materials, anodized aluminum, silicone and materials used in anesthesia utensils
- Good visibility of instruments in the solution
- Usable in all degrees of water hardness

Special properties:

- First-class cleaning performance
- Multi-enzymatic formula based on protease, lipase and amylase for an active residue-free cleaning action
- Free of boric acid and borates, free of perfumes and colorants

Application and dosage:

neodisher MultiZym can be used in immersion baths or ultrasonic baths. The dosing amount and application parameters depend on the area of application and the degree of soiling of the instruments. The following parameters are recommended when using neodisher MultiZym:

Application recommendations	
Manual cleaning / normal dirt load	0.125 oz/gal – 0.375 oz/gal, 59 - 122 °F, 2 - 10 min 1 – 2.5 ml/l, 15 - 50 °C, 2 - 10 min
Manual cleaning / high dirt load	0.375 oz/gal – 3.75 oz/gal, 59 - 122 °F, 2 - 10 min 2.5 - 30 ml/l, 15 - 50 °C, 2 - 10 min
Manual cleaning in ultrasonic baths	0.125 oz/gal – 3.75 oz/gal, 59 - 122 °F, 2 - 10 min 1 - 30 ml/l, 15 - 50 °C, 2 - 10 min

The neodisher MultiZym solution has to be rinsed off completely (preferably with deionized water).

The solution must be renewed at least daily; immediately in the case of visible soiling.

Notes on application:

- For professional use only.
- For economical and controlled dosing we recommend using manual dosing ancillaries, such as e.g. mixing equipment or dosing caps. Please contact us.
- Reprocessing should comply with all ordinances pursuant to the Medical Device Directive and should be performed with appropriate validated processes.
- Please observe the instrument manufacturer's instructions for use according to the requirements of public standards such as ANSI/AAMI ST 79:2017.
- Do not mix with other products.



neodisher® MultiZym



Technical data:

pH range	8.4 - 8.6 (0.125 oz/gal – 3.75 oz/gal 1 - 30 ml/l, determined in deionized water, 68 °F/20 °C)
Viscosity	< 50 mPa s (concentrate, 68 °F/ 20 °C)
Density	1.1 g/cm ³ (68 °F/20 °C)

Ingredients:


All components are listed in the TSCA inventory.
Isotridecanol, ethoxylated;
Concentration: ≥5<15%

Sodium alkylsulfonate;
Concentration: ≥1<5%

CE-mark:

neodisher MultiZym conforms to the European
directive 93/42/EEC, Annex I, concerning medical
devices.

Storage information:

Always store between 32 °F and 77 °F /
0 °C and 25 °C. Keep away from sunlight.
Usable for 2 years when stored as recommended.
For expiry date refer to the stamp mark on the
label behind the hourglass symbol .

Hazard and precautionary statements:

neodisher MultiZym is not a hazardous product
according to the Appendix D, OSHA Hazard
Communication Standard 29 CFR 1910:1200.
Dispose only when container is empty and closed.
For disposal of product residues, refer to Safety
Data Sheet.

MB 4087/3-5

Date of issue: 6/2019



DISTRIBUTED BY KEY SURGICAL LLC
8101 WALLACE ROAD / EDEN PRAIRIE, MN 55344 / USA
PH: 800.541.7995 KEYSURGICAL.COM

With the above information, which is appropriate to our current knowledge we describe our product regarding possible safety necessities, but we do not involve any quality description or promise certain properties.



Manufactured by:
Chemische Fabrik Dr. Weigert GmbH & Co. KG
Mühlenhagen 85 • Germany-20539 Hamburg

Phone: +49-40-7 89 60-0 info@drweigert.com
Fax: +49-40-7 89 60-120 www.drweigert.com