


FOAM CUBES AND SHEETS

Product Code(s)

- FOBxxx Foam Cubes
- FOSxxx Foam Sheets

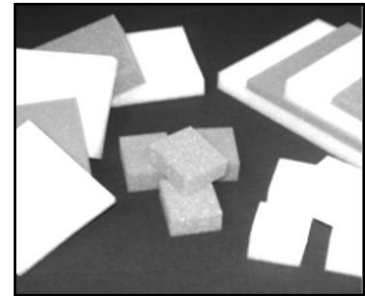

Foam Cubes and Sheets are supplied in any size and are unsterile. Available in standard white foam or charcoal grey, fine foam.

Material

Both standard white and fine charcoal grey foam cubes and sheets are manufactured from polyester-based polyurethane foams. They are specially formulated to be free of antioxidants that may be harmful in a medical environment.

Properties

The foam has an open cell structure which enables steam penetration through the cells during sterilisation. The foam is resistant to deformation during steam sterilisation.

Some discolouration may occur during steam sterilisation, but properties remain unaltered.

Specification

<u>Colour</u>	<u>White</u>	<u>Charcoal grey</u>
Density	24-30kg/m ³	20-22kg/m ³
Water Absorb.	35-40g/gram	35-40g/gram

Sterilisation

Foam Cubes and sheets may be sterilised by Steam Sterilisation up to 137°C/3.5min

Foam Cubes and sheets are not suitable for sterilisation by ethylene oxide or Plasma sterilisation techniques.

Storage Condition

Polyurethane foams should ideally be stored and used under moderate temperature and humidity conditions (ie.20°C, 50% relative humidity). Direct foam exposure to extremes of temperature and humidity can lead to accelerated foam degradation and discolouration.