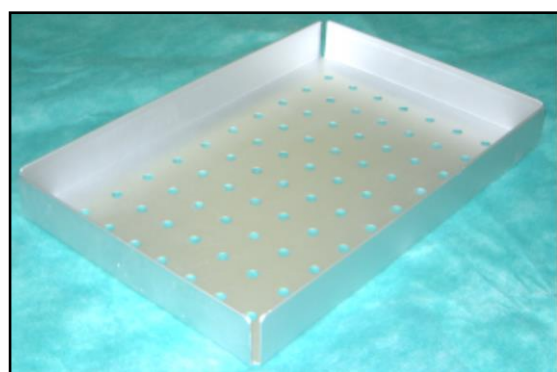
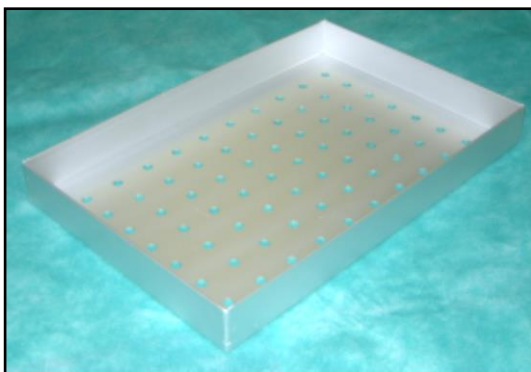


ALUMINIUM TRAYS**Product Code(s)**

➤ TRY001 – TRY999

A range of strong anodised aluminium trays suitable for steam sterilisation, constructed using 16 SWG aluminium sheets with folded edges for extra rigidity, they are provided with side holes for easy fixing of tray tags and are available with open or welded corners.

Lids are available for each size.

**Material**

Aluminium trays are manufactured from anodised aluminium 1050.

Cleaning and Sterilisation

Aluminium trays may be cleaned in a washer-disinfector and then sterilised in a steam steriliser at up to 137°C.

Ideally, an enzymatic cleaning concentrate should be used in the washer-disinfector. If this is not possible, the concentrate should be checked to ensure its compatibility with aluminium. Cleaning concentrates with high pH will discolour and corrode the aluminium.

Sizes

TRY001	-305 x 152 x 51mm
TRY002	-305 x 305 x 51mm
TRY003	-305 x 457 x 51mm
TRY004	-305 x 610 x 51mm
TRY005	-457 x 457 x 51mm with open corners