


**PLASTIC STERILISATION TRAYS****Product Code(s)**

➤ See Page 2

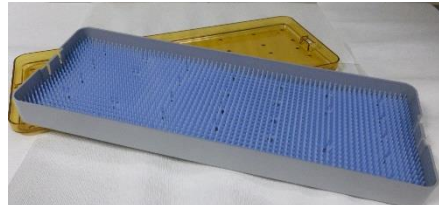

A complete range of plastic sterilisation trays for protecting instruments during sterilisation and subsequent transport. Special grid patterned holes provide drainage and fixing positions for a wide range of inserts. The tray bottoms have raised buttons to improve drainage during steam sterilisation. Single and double-decker trays are available.

**Materials**

Tray lids and bases are manufactured from a polyether imide engineering polymer. Silicone mats and bars are manufactured from medical grade silicone polymer.

**Designs**

Mini Trays  
Large microsurgical Instrument Trays  
Diamond Knife Tray  
Microsurgical Instrument Trays  
Scope Trays  
Deep Scope Trays

**Washing instructions**

Plastic sterilisation trays may be washed in any automated washer-decontaminator. The trays are compatible with all washer-disinfector cycles up to 95°C for a period of 20 minutes using any of the available detergents, including pH neutral detergents, currently on the market.

**Sterilisation**

Plastic sterilisation trays may be used in steam sterilisation up to 137°C/3.5 min, ethylene oxide and formaldehyde-based sterilisation systems.

**Plastic Sterilisation Trays**

Document Ref: CDS-3055

Page 1 of 2

Effective Date: 05 Oct 2019

TMP-CDS-V02: Effective 12 Feb 2018

Version No: 03

<b>CODE</b>	<b>Tray Type</b>	<b>Contents</b>	<b>Dimensions W x L x D (Inches)</b>
P0278	Mini Tray	Base, Lid & Mat	1.5 x 2.7 x 0.75
P0278125	Mini Tray	Base, Lid & Mat	1.5 x 2.7 x 1.25
P0625	Microsurgical Instrument Tray	Base, Lid & Mat	2.5 x 6.0 x 0.75
P062515	Microsurgical Instrument Tray	Base, Lid (no mat)	2.5 x 6.0 x 0.75
P0625K1	Diamond Knife Tray	1 Slot—Base, Lid, & 2 Bars	2.5 x 6.0 x 0.75
P0625K2	Diamond Knife Tray	2 Slots—Base, Lid, & 2 Bars	2.5 x 6.0 x 0.75
P0625K3	Diamond Knife Tray	3 Slots—Base, Lid, & 2 Bars	2.5 x 6.0 x 0.75
P2610A	Microsurgical Instrument Tray	Base, Lid & Mat	2.5 x 6.0 x 1.5
P0725	Microsurgical Instrument Tray	Base, Lid & Mat	2.5 x 7.5 x 0.75
P6540	Microsurgical Instrument Tray	Base, Lid & Mat	4.0 x 6.5 x 0.75
P7540	Microsurgical Instrument Tray	Base, Lid & Mat	4.0 x 7.5 x 0.75
P100675	Microsurgical Instrument Tray	Base, Lid & Mat	6.0 x 10 x 0.75
P100615	Microsurgical Instrument Tray	Base, Lid (no mat)	6.0 x 10 x 1.5
P1006151	Microsurgical Instrument Tray	Base, Lid, & Mat	6.0 x 10 x 1.5
P1006153	Microsurgical Instrument Tray	Base, Lid, Insert tray, 2 mats	6.0 x 10 x 1.5
P1500A	Large Microsurgical Instrument Tray	Base, Lid & Mat	10 x 15 x 0.75
P1015251	Large Microsurgical Instrument Tray	Base, Lid & Mat	10 x 15 x 1.5
P101525	Large Microsurgical Instrument Tray	Base, Lid, Insert tray, 2 mats	10 x 15 x 1.5
P261615	Scope Tray	Base, Lid, & Bars	2.6 x 16 x 1.5
P6520A	Scope Tray	Base, Lid, & Bars	3.0 x 22 x 1.5
P651815	Scope Tray	Base & Lid	6.5 x 18 x 1.5
P651815M	Scope Tray	Base, Lid, & Mat	6.5 x 18 x 1.5
P651830P3	Deep Scope Tray	Base, Lid, & Bars	6.5 x 18 x 3.0
P651830/HED	Deep Scope Tray	Base, Lid, & Mat	6.5 x 18 x 3.0