


## CLINIPAK MILLENIUM BLUE TRAY LINER

## Product Codes

- MILB001 – MILB999
- MILB000X000 – MILB999X999

Clinipak Millennium Blue Tray Liner is a non-woven material offering very high absorbency, making it ideal for use as a tray liner to absorb any excess moisture produced during sterilisation. Clinipak Millennium Blue is capable of absorbing 6-7 times its own weight in water, locking it in to its own porous structure preventing the moisture from adversely affecting the efficiency of the sterilisation process.

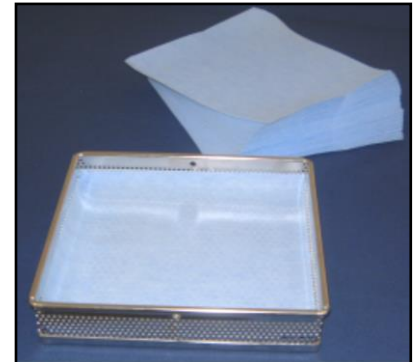

The material is supplied in a high visibility blue colour and possesses excellent mechanical properties to withstand handling and sterilisation.

The material is low linting and its very low weight to volume ratio makes it exceptionally easy to handle even in bulk.

## Sterilisation

Clinipak Millennium Blue Tray Liner is designed for use in steam sterilising systems at temperatures of 134°C and it is also compatible with ethylene oxide sterilising systems.

It is **not** suitable for use in plasma sterilising systems.

## Standard Properties

<u>Properties</u>	<u>Typical Value</u>		<u>Test</u>
<b>Weight</b>	40	Gsm	TM 1150
<b>Tensile Strength</b>			
MD	85	N/1 ply	TM 1178
CD	20	N/1 ply	TM 1179
CD wet	10	N/1 ply	TM 1187
<b>Thickness</b>	0,34	mm/1 ply	TM 1193
<b>Absorption Time</b>	<2	sec	TM 1015
<b>Absorption capacity</b>	750	%	TM 1016

The material is available in a range of sizes, optimised to fit most popular sizes of trays, though bespoke sizes are available on request.