


## SINGLE USE WASHER ID TAGS

### Product Code(s)

- TAG-WASH
- TAG-WASH-GRN, BLU, RED, ORA, YEL & PINK

Single-Use Washer ID Tags are designed for the identification of trays that are to be cleaned in a washer-disinfector.

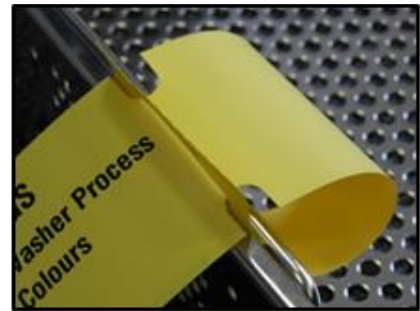

### Material

Single-Use Washer ID Tags are manufactured from high density polyethylene.

### Properties

Size	50 x 228.6 mm
Thickness	250 micron (except TAG-WASH which is 200 micron)
Weight	196 gsm
Colours	white, green, blue, red, orange, pink and yellow

### Printability

Variable information printing and barcodes can be put on both sides using thermal transfer printers. A black PEN001 may also be used to hand-write information on the tag.

### Washer-Disinfector compatibility

Single-Use Washer ID Tags are compatible with all washer-disinfectors that operate up to 95°C.

### Sterilisation Compatibility

Single-Use Washer ID Tags are not compatible with any sterilisation system and should be removed prior to any sterilisation procedure.