


CLINIPAK CHOICE BLUE ULTRA STERILISATION WRAP

Product Code(s)

➤ CBUxxxx (where xxxx= sheets size)

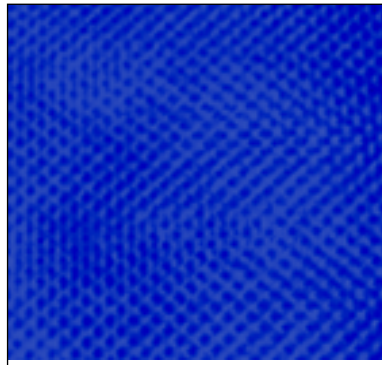

Introduction

Clinipak Choice Blue Ultra Wrap is a 70gsm, latex-free, 4-layer non-woven tray wrap, manufactured from spun-bonded and melt blown polypropylene, which provides a wet and dry microbiological barrier.

It is intended to be used as an outer sterilisation wrap where the very high strength provides protection against cuts and tears, particularly with heavy orthopaedic sets.

Clinipak Choice Blue Ultra Wrap is suitable for sterilization by Steam, Plasma (the STERRAD® system or a suitable alternative) and ETO. It can be used as an interleave in combination with other Clinipak barrier wraps – Choice and Choice Dry wraps as well as Plain and Crepe Papers.

Compliance

The material is compliant with EN 868 Part 2 and all relevant aspects of ISO 11607-1.

Typical Data

<i>Basis Weight</i>	<i>70 gsm</i>
<i>Tensile Strength MD</i>	<i>3.4 kN/m</i>
<i>Tensile Strength CD</i>	<i>2.1 kN/m</i>
<i>Elongation MD</i>	<i>70 %</i>
<i>Elongation CD</i>	<i>75 %</i>
<i>Microbiological barrier (wet)</i>	<i>Pass</i>
<i>Water Resistance (Mason Jar)</i>	<i>>75 min</i>
<i>Alcohol Repellency (70% alcohol)</i>	<i>Not alcohol repellent</i>

Storage/Shelf Life:

Clinipak Choice Blue Ultra Wrap has a shelf life of 5 years prior to use when stored at 20-80%RH and 10-35°C.