

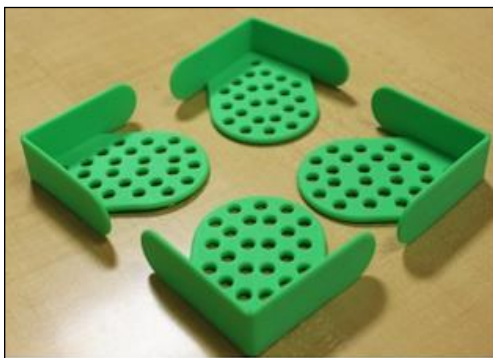
**SILICONE CORNER AND EDGE PROTECTORS****Product Code(s)**

- STOP-7100
- STOP-7100P500

Silicone Corner And Edge Protectors are designed to eliminate the risk of wrap damage caused by large, heavy sets or baskets with sharp corners. Manufactured from silicone, these Silicone Corner And Edge Protectors are multi-use. The corners are designed to be effective on trays with both straight and round corners.

With a large base containing holes, these Silicone Corners And Edge Protectors promote complete circulation of the sterilisation agent and also allow for better air circulation which in turn reduces drying times.

To use, simply slide them under the corners and/or feet of the basket prior to wrapping.

**Material**

Silicone Corner And Edge Protectors are manufactured from Silicone rubber.  
Latex Free.

**Washing and Sterilisation compatibility**

Silicone Corners And Edge Protectors are compatible with all types of washing disinfecting cycles.  
Silicone Corners And Edge Protectors are compatible with steam (up to 137°C), ethylene oxide and plasma sterilisation processes.

**Packaging**

STOP-7100                      40 individual corners per pack  
STOP-7100P500              500 individual corners per pack