

## CLINIPAK CHOICE GREEN X WRAP


## Product Code(s)

➤ CGX/xxxx (where xxxx= sheet size)

Clinipak Choice Green X Wrap is a latex-free, 4-layer non-woven tray wrap, manufactured from spun-bonded and melt-blown polypropylene, which provides a wet and dry microbiological barrier.

Clinipak Choice Green X Wrap is not alcohol repellent.

This material is suitable for sterilisation by Steam, Plasma (the STERRAD® system or a suitable alternative) and ETO and is intended for use as a wrap for trays & baskets. It can be used as an interleave in combination with other Clinipak barrier wraps – Choice, Choice Ultra and Choice Dry wraps as well as with Clinipak non-barrier Transportation Wraps and Plain and Crepe Papers.

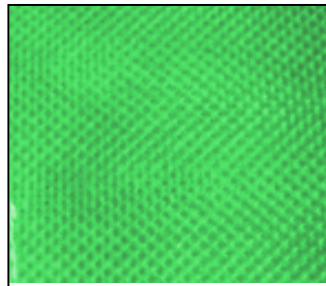

## Compliance

This material is in compliance with EN 868 Part 2 and all relevant aspects of ISO 11607-1.

## Standard Data

<i>Basis Weight</i>	55 gsm
<i>Tensile Strength MD</i>	2.0 kN/m
<i>Tensile Strength CDS</i>	1.1 kN/m
<i>Elongation MD</i>	50 %
<i>Elongation CD</i>	50 %
<i>Air Permeability</i>	220 l/min/100cm <sup>2</sup>
<i>Microbiological Barrier (dry)</i>	Pass
<i>Microbiological barrier (wet)</i>	Pass
<i>Water Resistance (Mason Jar)</i>	>75 min
<i>Alcohol Repellency (70% alcohol)</i>	Not Alcohol Repellent

## Sterilisation Compatibility

Clinipak Choice Green X Wrap is compatible with steam sterilisation processes up to 137°C, ethylene oxide sterilisation processes and plasma sterilisation processes.

## Storage/Shelf Life

Clinipak Choice Green X Wrap has a shelf life of 5 years prior to use when stored at 20-80%RH and 10-35°C.