


## SINGLE ENDED CLEANING BRUSHES

## Product Code(s)

Product Code	Sheath Diameter	Brush Length	Diameter	Head length	Body
BBRW60SE5	3.3mm	60cm	5.5mm	45mm	Twisted Wire
BBRW60SE8	3.3mm	60cm	8mm	45mm	Twisted Wire
BBRP60SE5	3.3mm	60cm	5.5mm	45mm	Purple Sheath
BBRP60SE8	3.3mm	60cm	8mm	45mm	Purple Sheath
BBRP60SE11	3.3mm	60cm	11mm	45mm	Purple Sheath
BBRP60SE14	3.3mm	60cm	14mm	45mm	Purple Sheath
BBRP60SE18	3.3mm	60cm	18mm	45mm	Purple Sheath

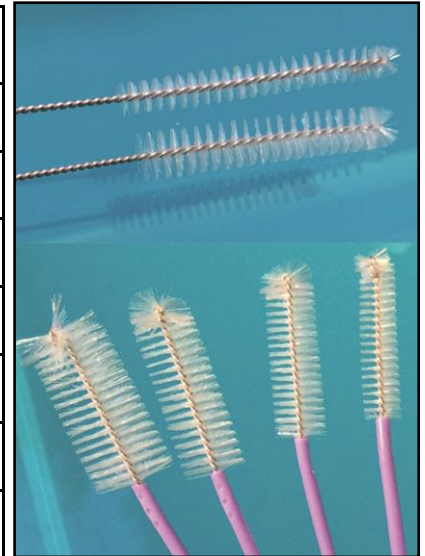

A range of single-use brushes for use in association with water and other solvents to remove organic soils and other residuals from channels and valves of channelled rigid instruments and scopes. Brushes are latex-free.

**Suitable for Cleaning Rigid Instruments/Scopes**

Laparoscopic and Orthopaedic, Suction Tubes and Guides Cystoscopes, Resectoscopes and Arthroscopes

**Material & Properties:**

BBRW60SE(X) – Nylon head with twisted wire body

BBRP60SE(XX) – Nylon head with purple Delrin® sheath body

Delrin® is more rigid than polypropylene so it is easier to insert and manipulate. Delrin® based products are used in precision parts which require high stiffness, low friction and excellent dimensional stability. These strengths in turn reduce the time and labour required to clean items.

**High Density Bristles**

Bristle fan tip protects rigid instruments/scopes; no sharp ends or hard metal tips that can scratch or cause damage. Available in a range of cleaning head sizes.

**Cleaning the Brushes**

BBRW60SE(X) and BBRP60SE(XX) Brushes are single use only.

BBRW60SE(X) and BBRP60SE(XX) Brushes cannot be re-used or re-processed.

*Note: Cleaning with the aide of brushes should be carried out in accordance with instrument /scope manufacturers' recommendations.*