




# T-EZ PRO™ - CLEANING INDICATOR


The T-EZ Pro contains a sheep blood compound that correlates to human blood. It is easy to use: just remove the T-EZ Pro from its package and attach it directly to the sieve tray. Thanks to the transparency of the plastic, there is no need to remove the indicator strip from its PCD for interpretation. However, it is possible to open the test body and archive the indicator for documentation purposes.

## ABOUT:

- Easy to use
- All in one
- Possible to remove of the indicator from its PCD for protein tests
- Single use
- Sheep blood compound as a verification indicator

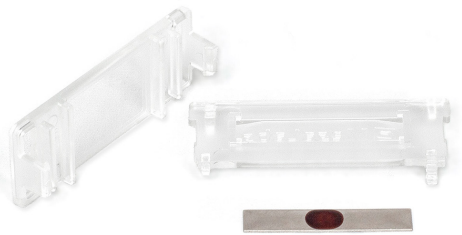

REF	DESCRIPTION	QTY
EZ-015-100	T-EZ Pro Cleaning indicator to test the cleaning performance of WD, tested according to ISO 15883	15 pcs

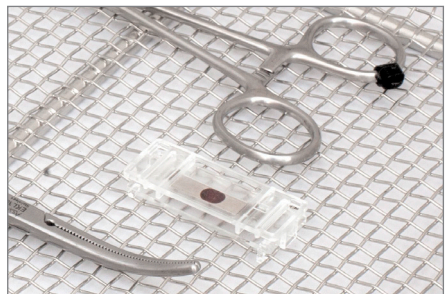

Manufactured by  
**MesaLabs**