


CHIROKLIP

Wire and Pin Guard



ChiroKlip™ is a protective guard used to protect the ends of exposed wires and pins to minimize the risk of 'needle stick' injuries and to aid in the removal of K-wires. A one-size-fits-all solution and alternative to traditional pin covers.

- Intended use with efficiency in mind, ChiroKlip fits most commercially available K-wire and pin diameter sizes (typical diameters range from 0.5mm to 4mm)
- Small, robust aluminum clip (butterfly-shaped) with strong adhesive inside that can fold around the exposed/protruding end of a K-wire
- Prevents 'needle stick' injuries for patients and staff
- Aids in the removal process of K-wires
- Designed to be disposed with the K-wire
- Sterile, ready-to-use
- Single use

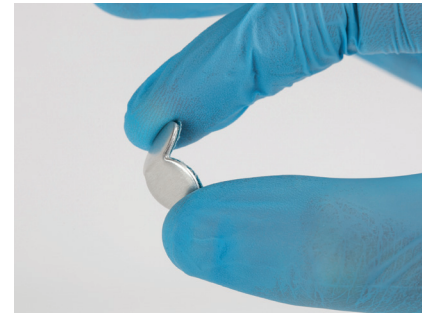

CLINICAL STUDY



Product No.	Description	Quantity
CK100	ChiroKlip Wire and Pin Guard	25/box (2 per sterile pack)

VISIT [KEYSURGICAL.COM](https://www.keysurgical.com) OR CALL 800.541.7995 FOR MORE INFORMATION