




ENDOSCOPE MANUAL CLEANING: GEMINI BRUSH™

The Gemini Brush is a double-ended cleaning brush designed for removing soil and debris from flexible endoscope channels during the manual cleaning process. Other features include a non-slip sheath which helps with grip during use as well as protective tips at the end of the brush head to protect the instrument channels from inadvertent damage. The Gemini Brush heads are tapered to support your need to clean a wide range of endoscopes.

FOR CHANNEL
SIZES
2.0mm - 7.0mm

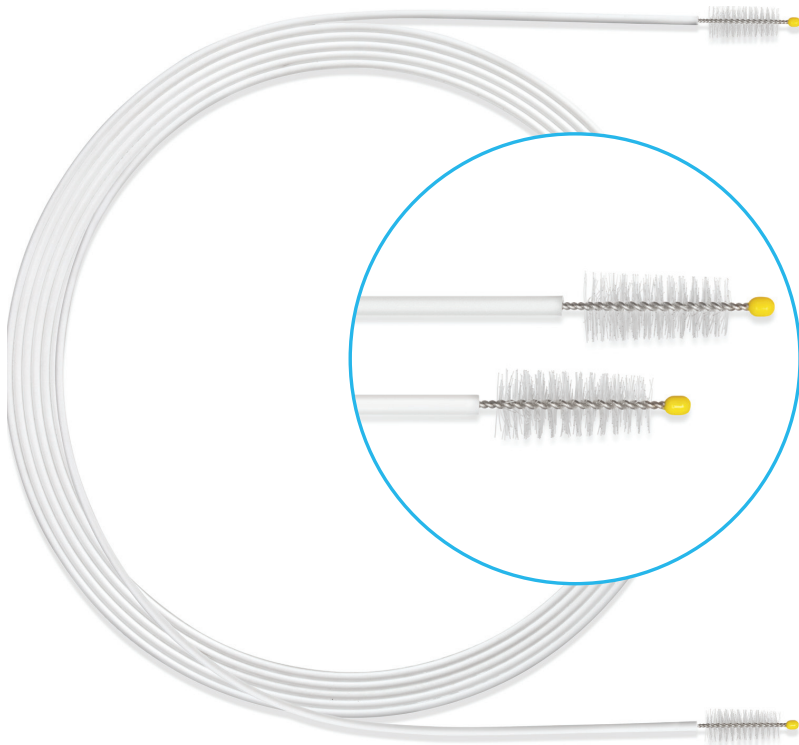

ABOUT:

- Gemini Brush features two tapered brush heads (one on each end) to push soil through the channel with the leading brush head and pull remaining soil through with the other
- Tapered brush heads (5-7mm) for cleaning scope channel diameters of 2.0mm up to 7.0mm
- Overall length: 230cm
- Textured, non-slip sheath for maintaining grip during use
- Protective tips help prevent inadvertent damage to the scope channel

PRODUCT NO.	DESCRIPTION	QUANTITY
BX00711607	Gemini Brush, Double-Header Channel Cleaning Brush, Tapered Brush Heads 5-7mm	50
BX00711627	Gemini Brush, Double-Header Channel Cleaning Brush, Tapered Brush Heads 5-7mm	100

VISIT [KEYSURGICAL.COM](https://www.keysurgical.com) OR CALL 800.541.7995 FOR MORE INFORMATION