



NEW! T-EZ PRO™: WASHING INDICATOR

Testing your washer(s) not only at first install but during daily use is vital to the entire process of surgical instrument cleaning performance. To challenge the function of your washer, use T-EZ Pro to test the effective removal of surgical soil from instruments and medical devices during the mechanical cleaning process. T-EZ Pro uses a standardized true blood test soil as a verification indicator, an accurate representation of the most common contamination found in healthcare.

“Mechanical cleaning equipment performance should be tested each day it is used and all results should be recorded.”


- AAMI ST79 7.6.4.5

“Employing a test device that is a consistent and repeatable challenge to the cleaning effectiveness of the equipment is critical to a successful quality assurance program.”

- AAMI ST79 13:2

ABOUT:

- T-EZ Pro is a single-use PCD manufactured in accordance with ISO/TS 15883-5:2005. Each PCD includes a clear plastic case and a stainless steel plate with test soil
- Uses true blood as a verification indicator
- Easy to view test plate securely placed in plastic holder
- Easy to see pass/fail results - with accurate calibration of time, temperature, water pressure, detergent concentration, and cycle parameters the test soil will be removed completely at the end of the cleaning cycle
- Supports AAMI ST79 best practices, dedication to patient safety, and overall workflow
- Each T-EZ Pro is packaged in a sealed aluminum pouch indicating product, lot number, and expiration date

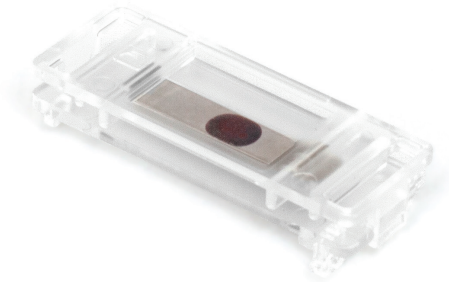

EASY TO READ PASS/FAIL CHART

PASS

FAIL



(NO REMAINING TEST SOIL)



PRODUCT NO.	DESCRIPTION	QUANTITY
EZ-015-100	T-EZ Pro cleaning indicator test	15/box



a STERIS company

VISIT OUR WEBSITE KEYSURGICAL.COM OR CALL 800.541.7995 FOR MORE INFORMATION

